


# PLANER

## DeltaTv7

### NEW SOFTWARE FOR YOUR CONTROLLED RATE FREEZER

- Stay informed with email alerts and real time profile updates during the freezing run
- Control multiple freezers with one instance of DeltaTv7 (Batch Mode)
- Remove the requirement for paper print outs with a comprehensive digital report for each freezing run
- Access your freezing run remotely

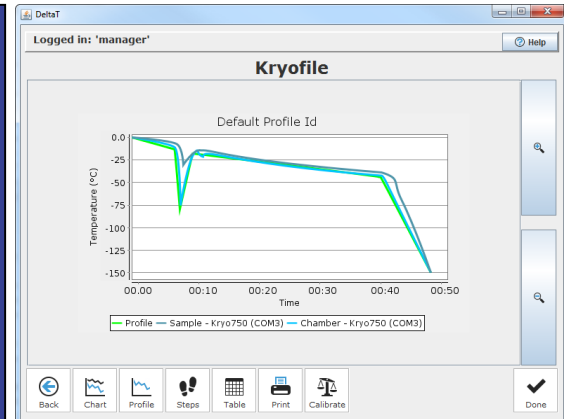

#### MULTIPLE FREEZER CONTROL

DeltaTv7 allows users to deploy one profile to several freezers. This means large numbers of samples can be run on different machines with the subsequent freezing profile for all samples being displayed on the same graph.

#### COMPLIANCE

Every freezing run produces a detailed date and time-stamped summary, which can automatically be stored to a secure location over a local area network. The inclusion of password protection and tamper-proof files permits scientists to use DeltaTv7 as part of a 21 CFR Part 11 solution.

#### CONTEXT

Users of DeltaTv7 can include additional context with the freezing profile with the ability to add time-stamped annotations both during and after the freezing run.

#### COMPREHENSIVE CONTROL AND VISIBILITY

Profiles can be created with up to 32 different stages (heat/cool/hold) combined with real-time updates of both the sample and chamber temperatures. Users are also provided with email notifications during the freezing run (in addition to the audible alarm of the freezer) in the case of significant deviations from the selected profile.

**For more details, please visit:**

**<http://www.planer.com/products/cryo-freezers/software/deltat-v7.html>**